



## SIBIONICS GS1 CGM System User Guide

[Home](#) » [SIBIONICS](#) » [SIBIONICS GS1 CGM System User Guide](#) 


# SIBIONICS

GS1 CGM System  
Quick Start Guide

### Contents

- [1 Get Started](#)
- [2 Set Up the App](#)
- [3 Get Ready](#)
- [4 Use SIBIONICS App](#)
- [5 Replace your Sensor](#)
- [6 Documents / Resources](#)
  - [6.1 References](#)
- [7 Related Posts](#)

### Get Started

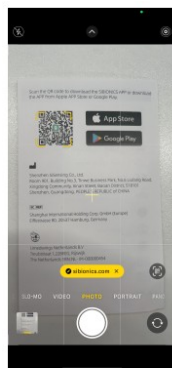



<https://guide.sibionics.com/cgm-app/>

Please use your phone camera to scan the QR code above to download the SIBIONICS App.

## Set Up the App



### Step 1:

Sign up with your registered email or follow the instructions to register an account with your email.

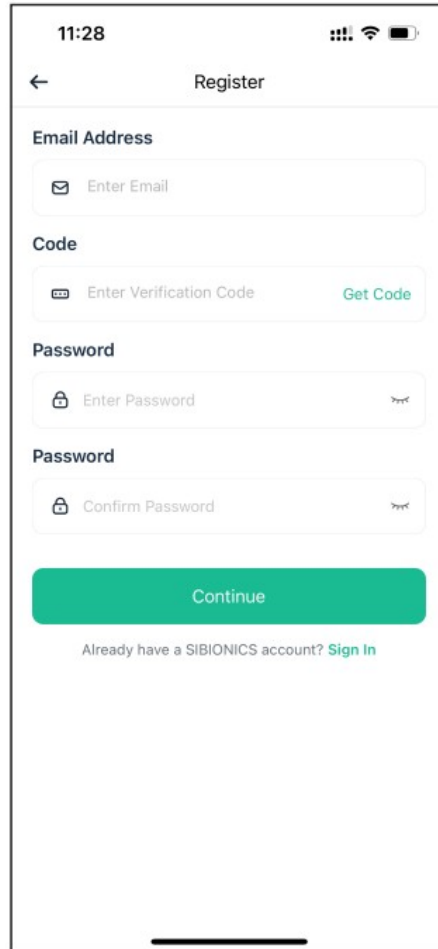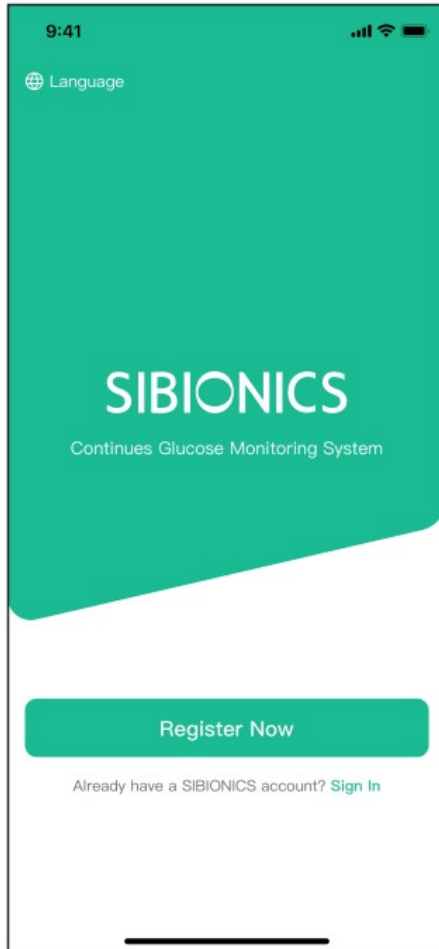


### Step 2:

Follow the App instructions to set up your App.

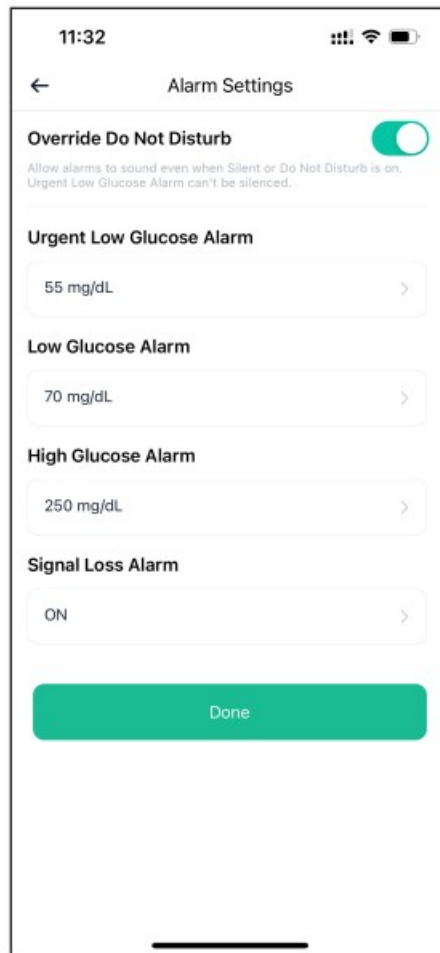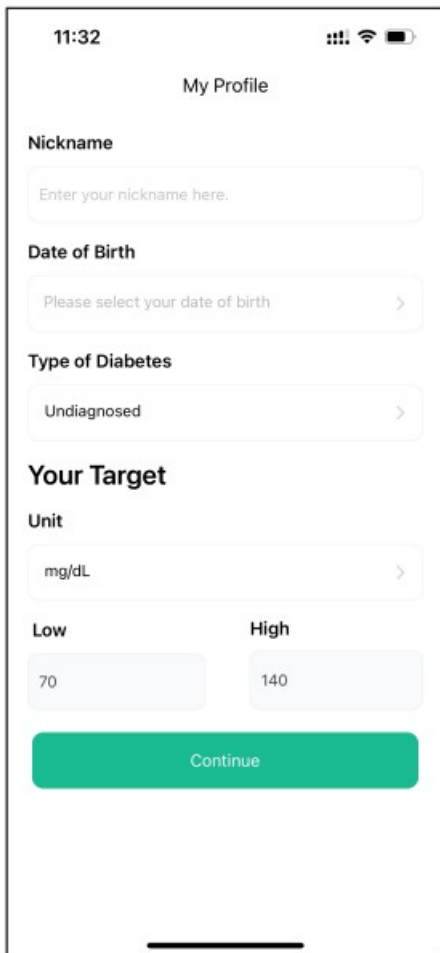


## Get Ready

**Step 1:** Select a site on the back of your upper arm.

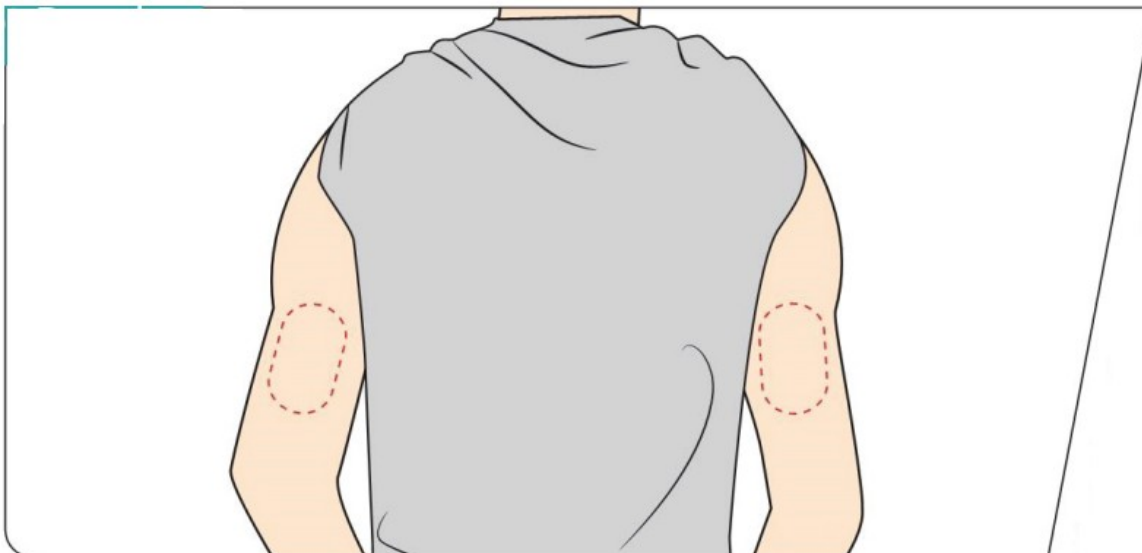


**Step 2:** Use an alcohol wipe to remove any oily residue and allow your skin to dry before applying the sensor.

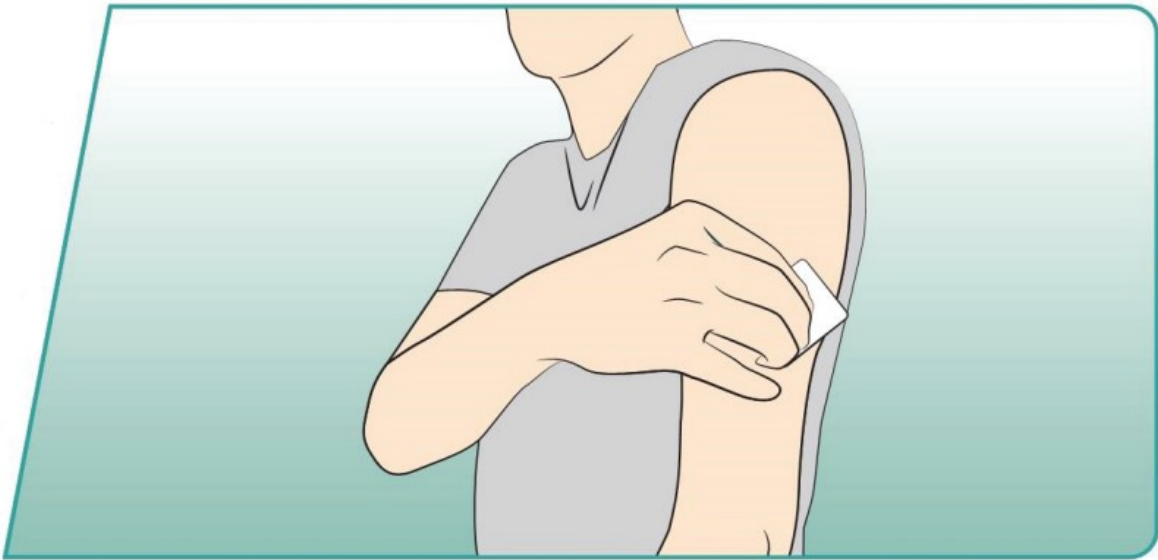


**Step 3:** Peel off the lid on the sensor pack.



**Caution:** The pack should be put on a flat desktop.

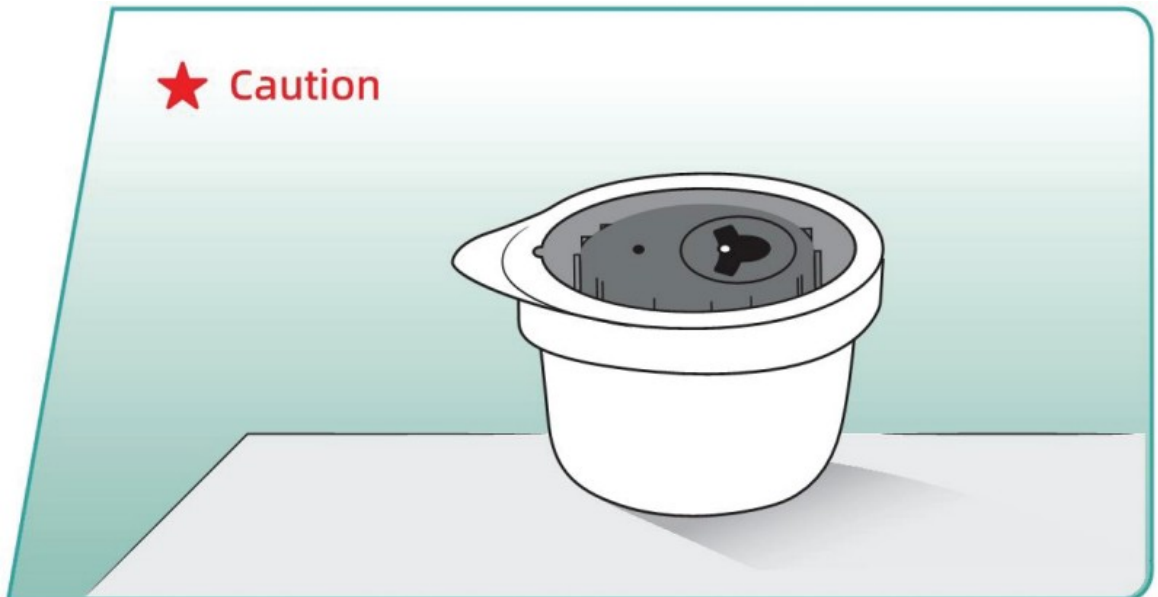


**Step 4:** Uncover the cap from the sensor applicator.

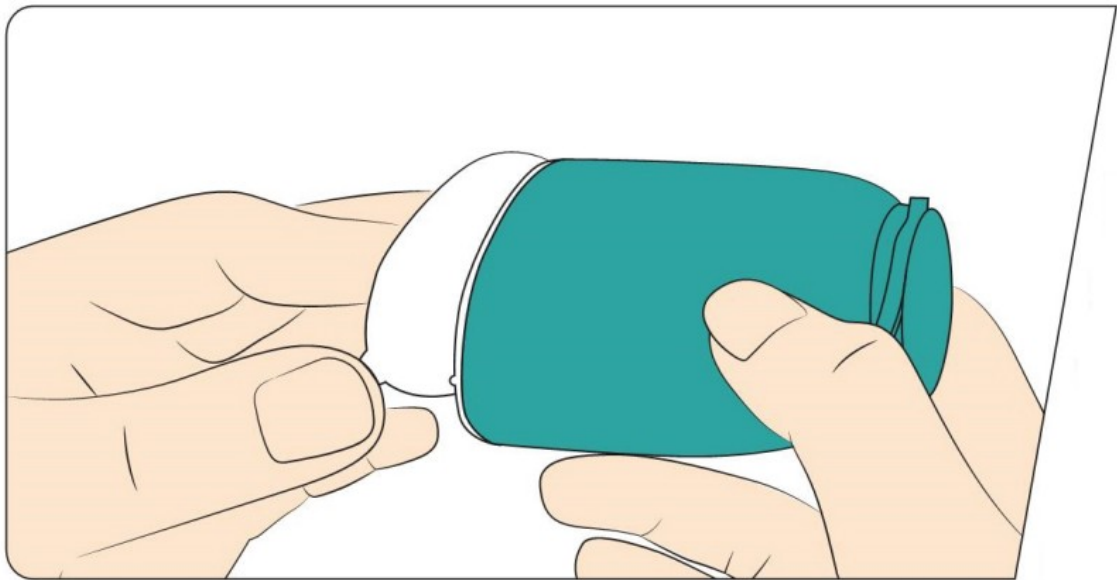


Press the sensor applicator firmly down on the pack until it comes to a stop.




**Step 5:** Lift the sensor applicator out of the sensor pack.

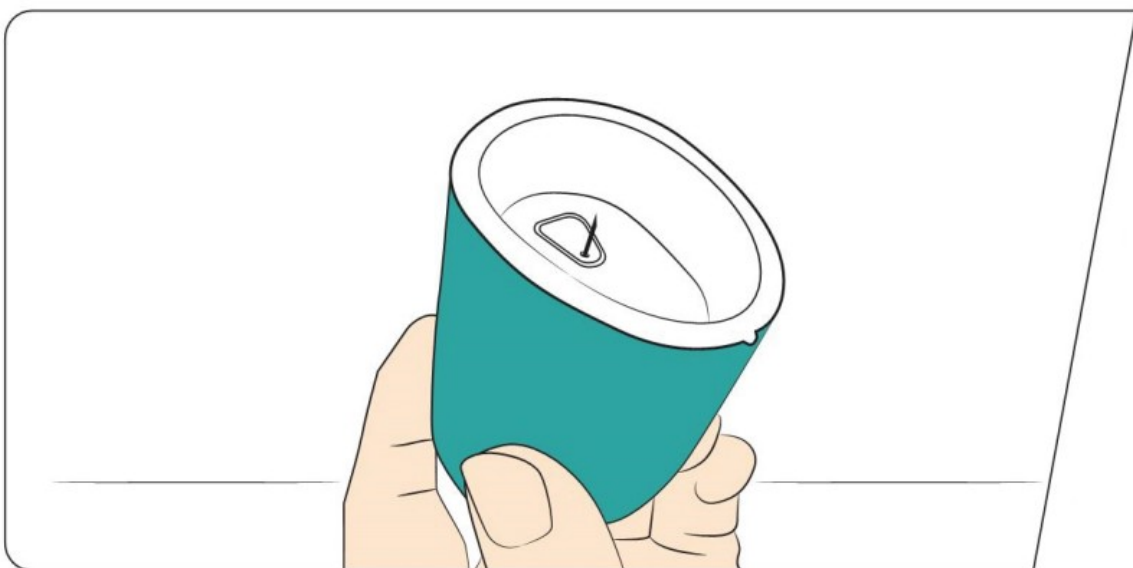


If the sensor is not ready to apply, please repeat Step 4.

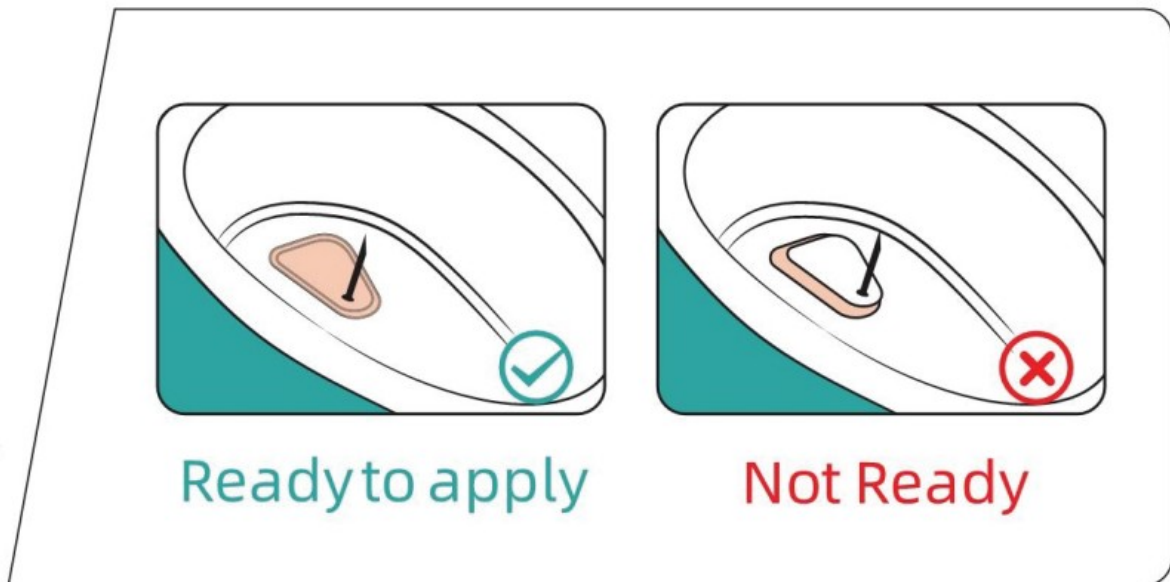


**Step 6:** Press and pull the safety clip out.

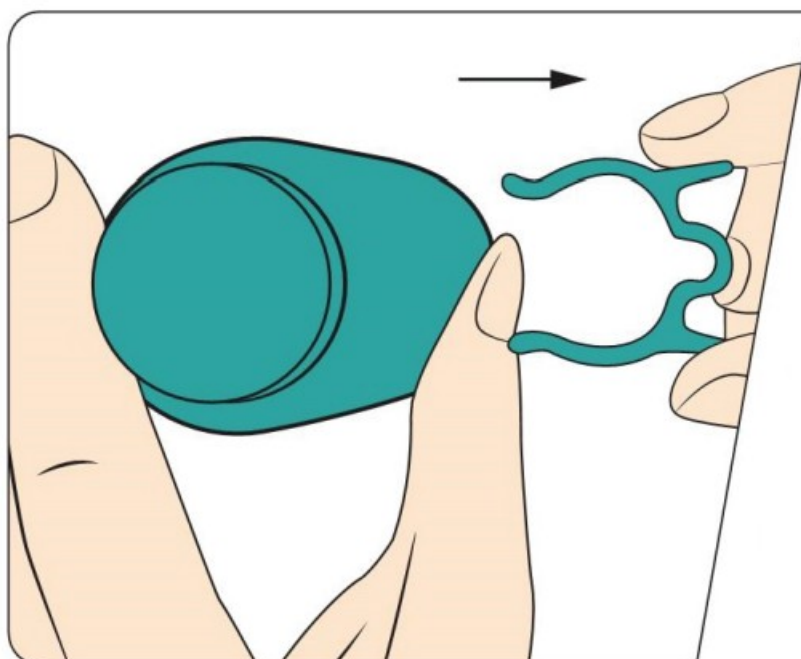

**Step 7:** Place the sensor applicator over the site and push the button down.

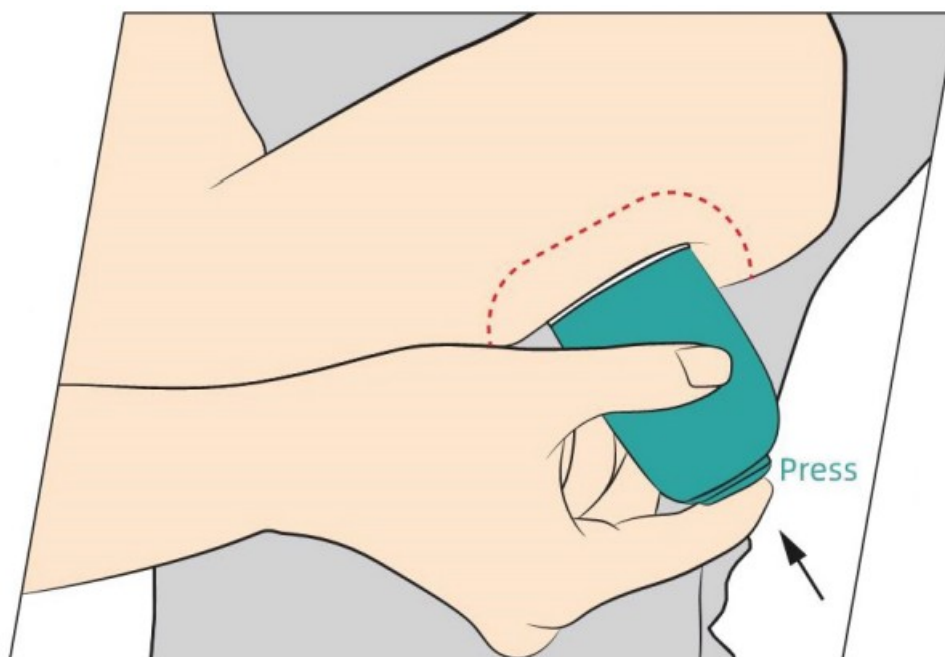


Now your sensor is well applied.

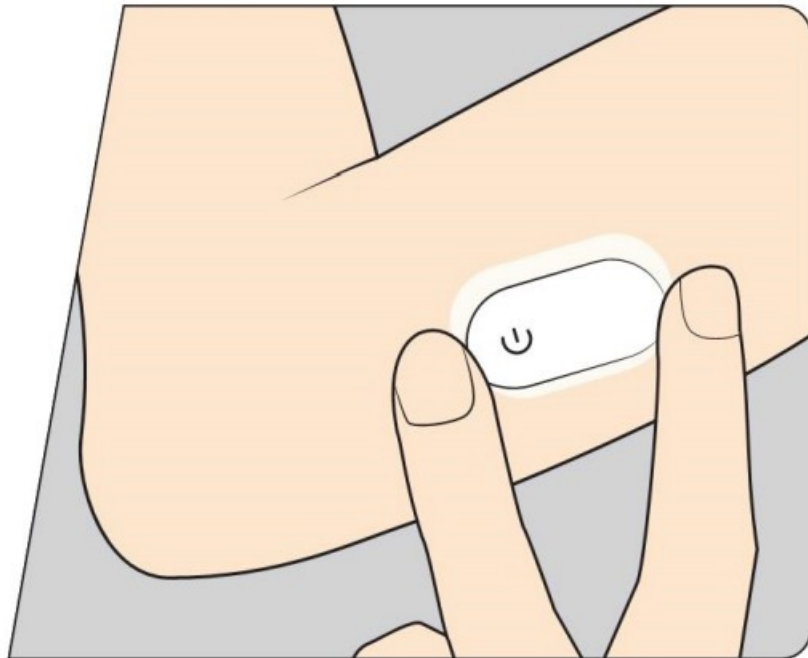


## Use SIBIONICS App

Open the App and click the button “Scan the code to connect”.


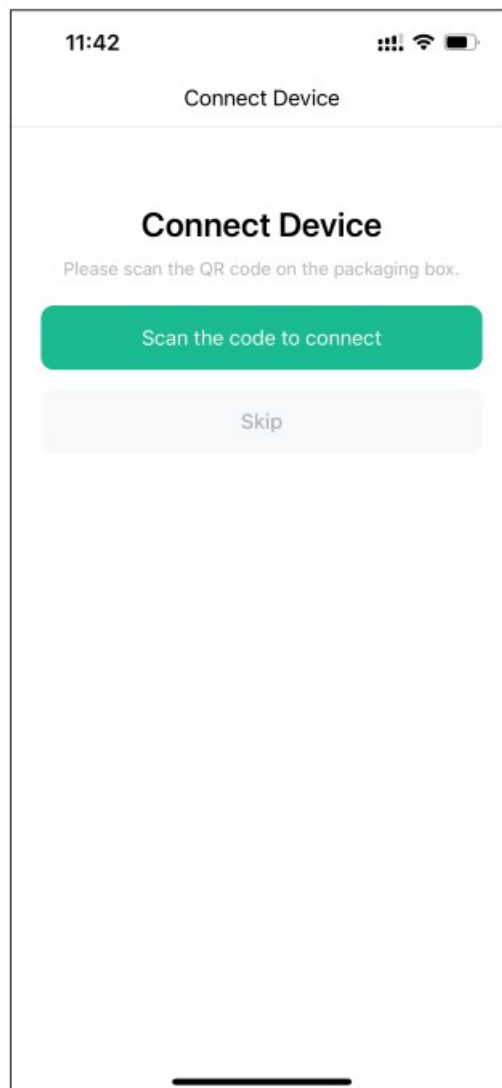


Use the camera to scan the QR code on the box.

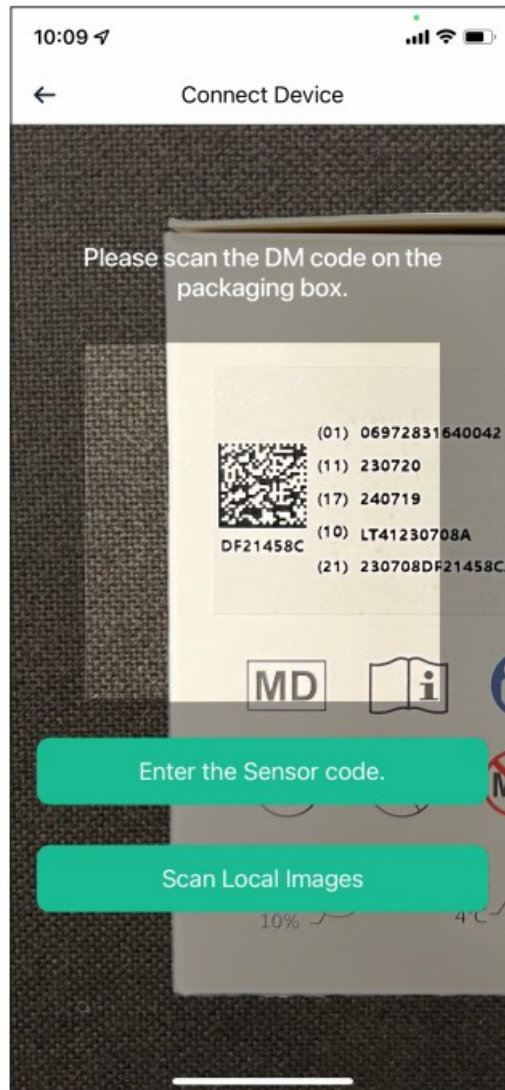

Sensor warm up timer tells you when you will start getting readings.



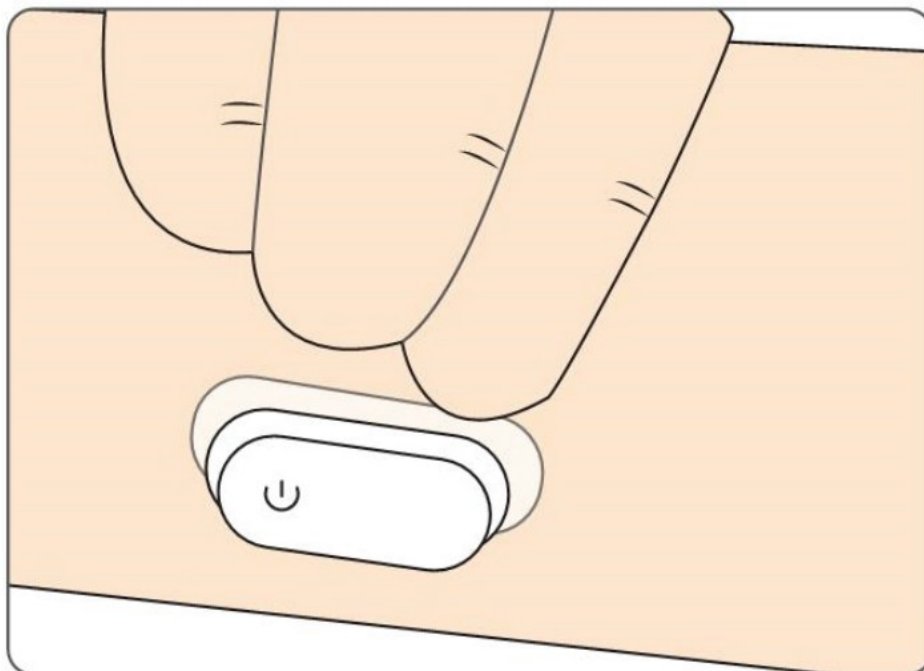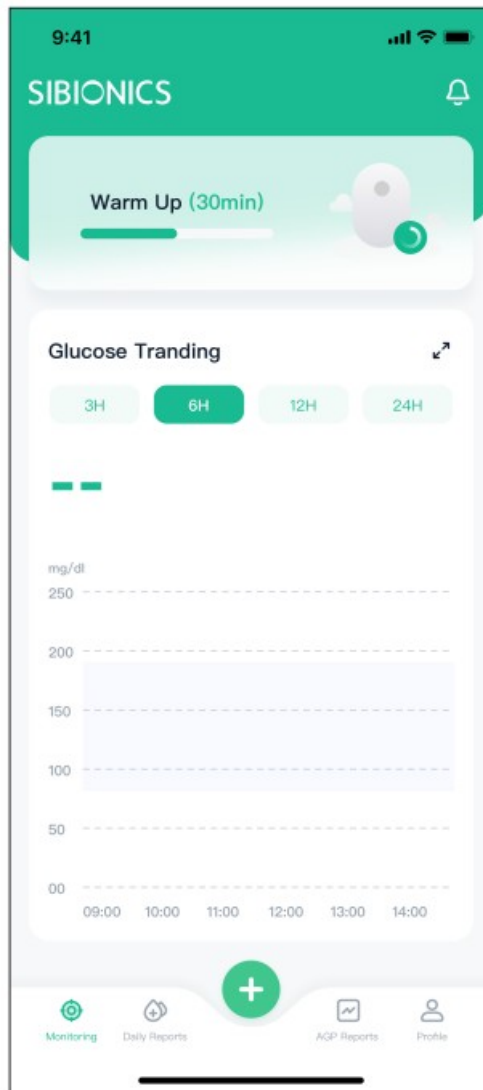

## Replace your Sensor


Your sensor automatically stops working after 14 days and must be replaced. You should also replace your sensor if you notice any irritation or discomfort at the application site.

Now you complete the quick guide and can explore the app by yourself now. Go to the User Guide to find out more.

# SIBIONICS

---

## Documents / Resources

	<p><a href="#">SIBIONICS GS1 CGM System</a> [pdf] User Guide GS1 CGM System, GS1, CGM System, System</p>
---	--

## References

- [User Manual](#)